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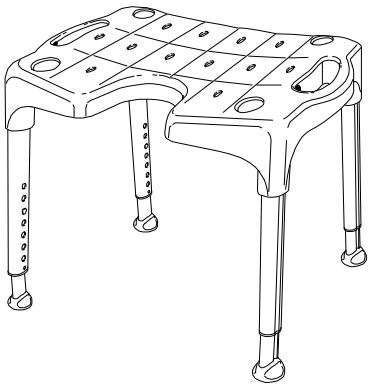
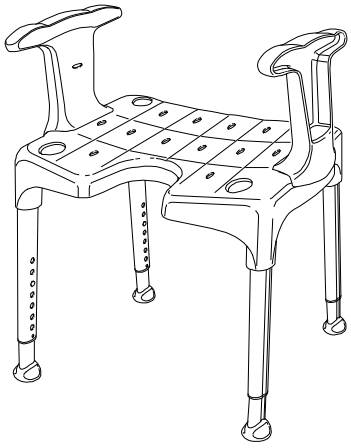
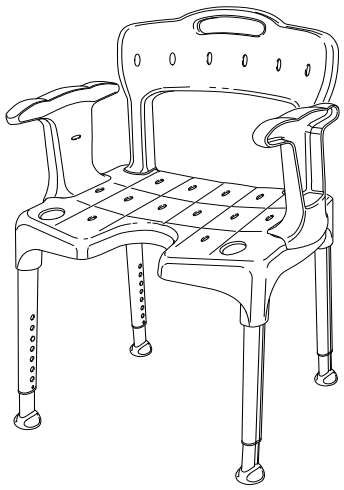
English

Thank you for choosing an Etac product. In order to avoid damage during assembly, handling and use, it is important to read this manual and save it for future reference. You can also find it at www.etac.com. You can select your language via the "International" and "Local websites" link. Here you will also find other product documentation, such as prescriber information, pre purchase guide, and reconditioning instructions. In the manual, the User is the person sitting in the chair. The carer is the person helping the user.

This symbol appears in the manual along with text. It indicates instances where the safety of the user or the assistant may be at risk.

The images which are referred to on each page are found at the start of the manual (pages 2 to 3). See Figures A to F for the corresponding text.

Description

Etac shower stools/shower chairs:

Etac's shower stools/shower chairs have a stable design with built-in flexibility that adapts to uneven floors. All products have height-adjustable aluminium legs that are simple to attach. The legs are angled outwards to maximise stability regardless of height. All products are equipped with ferrules that provide good grip, even on uneven and damp floors. The products can be tilted forward to further facilitate standing up.

The seats are designed with an anti-slip pattern, ensuring a safe and non-slip seat even when it is wet. The seats also have drainage holes that effectively remove excess water from the stool.

Swift shower stool/shower chair:

Etac Swift is available as a shower stool, shower stool with arm supports, and complete shower chair with arm supports and backrest. Arm supports/Backrest are fitted and removed without tools and can be added afterwards/removed if necessary. The seat is reversible and features a hygiene recess. If the arm supports are not used, the opening enables a stable grip and support.

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This product should be recycled as per national regulations.

Product description — See Figure A

	Swift	Smart	Easy	Edge
A1. Seat **A2.** Legs **A3.** Ferrule	●	●	●	●
A4. Product label*	●	●	●	●
A5. Armrest **A6.** Backrest	●			

The product's date of manufacture can be read from the barcode on the product. The number 11 is shown below the barcode in brackets. The number combination after these brackets is the date of manufacture.

Technical data — See Figure B

Assembly — See Figure C

Dismantling — See Figure D

1. Backrest 2. Legs 3. Armrest

Guarantee, lifetime

5-year guarantee against material and manufacturing defects. For terms and conditions, see www.etac.com
Service life 10 years. For complete information regarding the service life of the product, see www.etac.com.

Safety & Maintenance — See Figure E

- A defective product or accessory must not be used.

E1. A height difference of more than one hole (3°) is not allowed for chair legs.
E2. When moving sideways, do not lean on just one armrest.
E3. Max user weight
E4. Clean the product using a solvent-free cleaning agent with a pH value of 5-9 or with a 70% disinfectant solution.
E5. Can be decontaminated at max. 85° C/185°F for 3 minutes if required.

In case of an adverse event occurred in relation to use of the device, it should be reported to your local dealer and the national competent authority in a timely manner. The local dealer will forward information to the manufacturer.

Reuse

The product is suitable for reuse. The number of times it can be reused depends on how often and in which way the product is used. Before reuse, thoroughly clean and disinfect the product and check its safety features. If necessary, take remedial action.

Storage

The product should be stored indoors in a dry place at a temperature above 5°C. If the product has been stored for a long time (more than four months), its function must be checked by an expert before use.

Accessories — See Figure F

		Swift	Smart	Easy	Edge
F1.	Back support	●			
F2.	Armrest/pair	●			
F3.	Seat pad	●			
F4.	Backrest pad	●			
F5.	Soap dish/Shower head holder	●			
F6.	Seat pad		●		
F7.	Seat pad			●	
F8.	Seat pad				●
F9.	Swivel plate		●	●	

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Etac Supply Center AB