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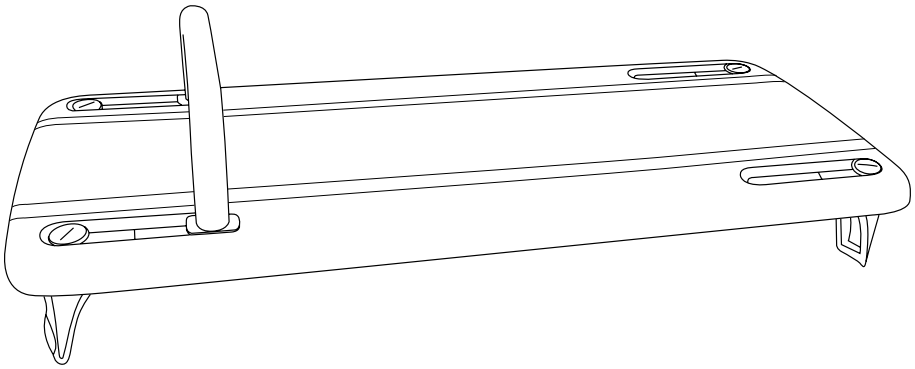
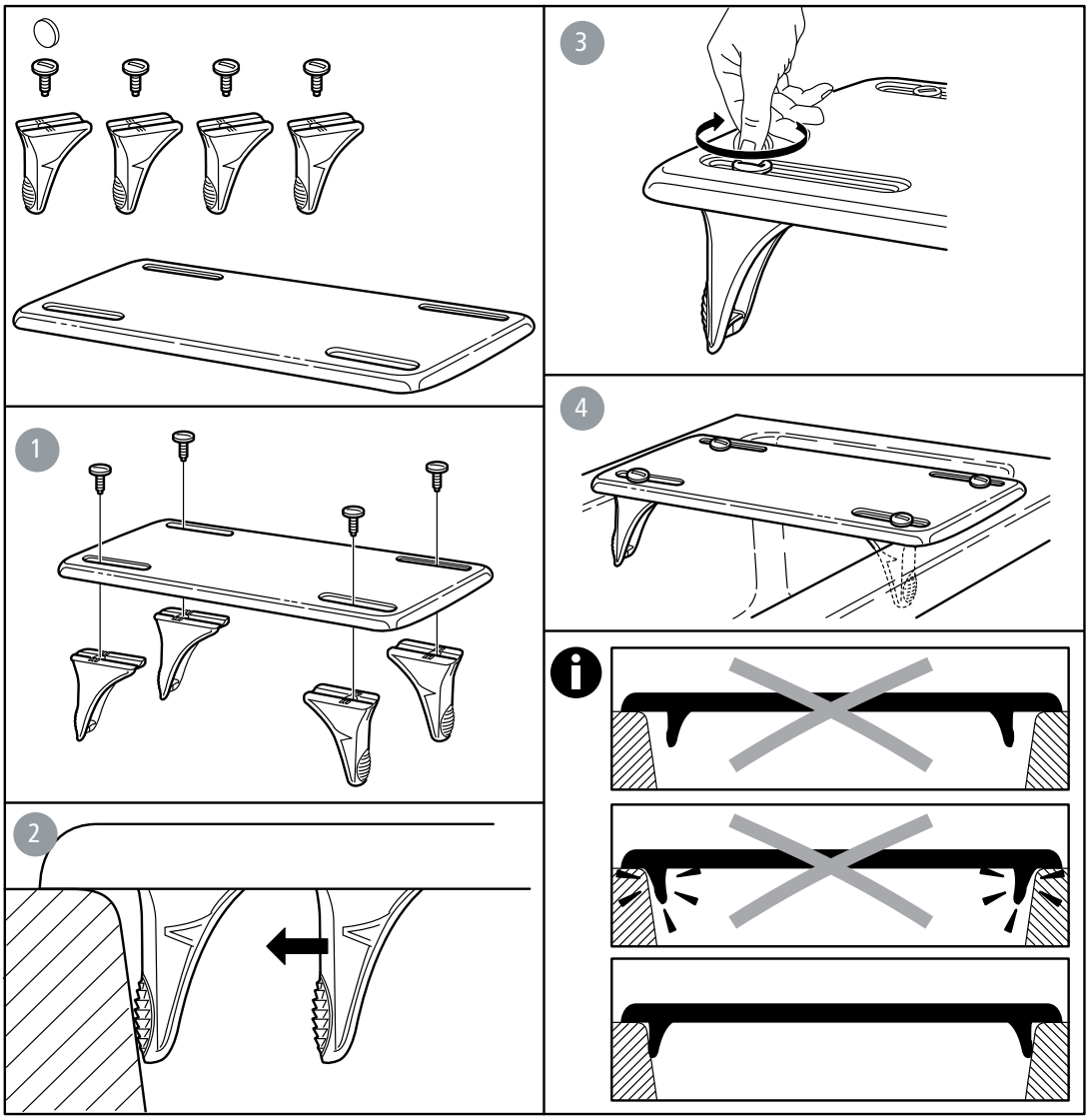
en English

Intended use

The Etac Rufus Plus bath board is designed for use as a bathtub seat. The bath board facilitates movement to the bathtub by offering a safe seated transfer over the side of the bath. The Etac Rufus Plus bath board is designed for persons with restricted mobility and a maximum user weight of 200 kg.

The product conforms to the requirements of the Medical Devices Regulation (EU) 2017/745.

Installation, use ***Figure A***
Maintenance and safety ***Figure B***

Max user weight: 200 kg.

Cleaning: Clean the product using a solvent-free cleaning agent with a pH value of 5-9 or with a 70% disinfectant solution.

Decontamination: Can be decontaminated at max. 85° C/185 °F for 3 minutes if required.

The materials used are corrosion-resistant.

The materials used are resistant to all general disinfectants.

Warning: A defective product must not be used.

In case of an adverse event occurred in relation to use of the device, it should be reported to your local dealer and the national competent authority in a timely manner. The local dealer will forward information to the manufacturer.

Guarantee, service life

5-year guarantee against material and manufacturing defects. For terms and conditions, see www.etac.com.

Service life 10 years. For complete information regarding the service life of the product, see www.etac.com.

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Etac Supply Center AB
Etac R82 GmbH
Etac R82 B.V.
Etac Ltd.