FLIPPER

Art. 271740, 271744, 271748, 271751, 271940, 271944, 271948 & 271951

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Mobilex A/S
Grønlandsvej 5
DK-8660 Skanderborg
www.mobilex.dk

1. Introduction

Thank you for buying a high-quality Mobilex product. These instructions for use contain a description and important information about the use of the medical device. It is intended to help you to use the product correctly and safely. It is very important to read the instructions carefully before using the device. Pay particular attention to the safety guidelines and always follow them. Mobilex is constantly trying to improve and upgrade its products. Therefore we reserve the right to change specifications and technical data of the product without further notice.

Intended use of the device

Intended use of the wheelchair is to move around in a sitting position, independently or with the help of an attendant. Using the wheelchair for purposes other than those indicated above, is prohibited.

Indications

Manual wheelchairs are indicated for people with mobility limitations who have problems with walking, standing or sitting without stable support. The wheelchair is indicated for use to move around in a sitting position, independently or with the help of an attendant. Using the wheelchair for purposes other than those indicated above, is prohibited.

Contraindications

The wheelchair should not be used by people who cannot maintain a stable sitting position. Contraindications cannot be related to one specific clinical picture, and they depend on individual evaluation of potential user's health status. The choice of orthopedic aid appropriate for a given clinical picture must be always made with the help of a specialist physician or physiotherapist. In addition, observe the safety information in Chapter 2.

Areas of use

The product can be used in home care, health and long-term care facilities, and hospitals.

Medical incidents and unwanted events

In case of any unwanted events related to the device, please contact the dealer or manufacturer:

Mobilex A/S		
Grønlandsvej 5	Phone:	+45 87 93 22 20
DK-8660 Skanderborg	Telefax:	+45 87 93 17 77
Denmark	Email:	info@mobilex.dk

2. Safety requirements

- Make sure that this manual has been read by all persons using the device.
- The manufacturer assumes no responsibility for damage and/or injury resulting from failure to follow the operating instructions.
- Only use the product in flawless condition.
- If defects or errors are detected you must immediately contact the dealer.
- Follow all manufacturer's recommendations and warnings.
- Be aware of the information on the product label.
- Only use the product for the described purpose.
- Don't climb on to the device
- Never make changes to the wheelchair.
- The device must only be used on stable surface.
- When adjusting the device the user must take care not to squeeze the fingers or other body parts.
- Ensure that all 4 wheels are always in contact with the surface during use. This will ensure proper wheelchair balance and will help avoid accidents.
- Do not exceed the maximum weight capacity of 136 kg.
- Using the wheelchair or exposing it to direct and prolonged contact with the sun may lead to skin burns.
- It is not recommended to use the wheelchair on the beach, in mud or in extreme weather conditions.
- The person assisting the disabled person in the wheelchair should be physically capable and completely healthy.
- Avoid leaving the wheelchair near sources of heat and open flames (campfires, stoves). Smoking is also not recommended during use as the wheelchair is not fire resistant.

Parts supplied

The device is supplied in a carton box, disassembled for easier transportation.
The contents of a package are:

- ✔ 1 main frame with installed: upholstered back and seat, 2 armrests, 2 drive wheels, 2 brakes, 2 steering wheels, 2 quick release axles, 2 anti-tips (optional);
- ✔ 2 detachable footrests;
- ✔ 1 user manual;
- ✔ 1 seat cushion (optional)


Preliminary check
Check the contents of the box with the list above. If at this stage you discover that a part is missing or damaged, contact your supplier immediately.

3. Description

Wheelchair provides:

- A comfortable transportation and storage thanks to its crossed frame structure and wheels with a quick-release axle
- Easiness of maneuvering
- Easy transfer to - and out of the wheelchair
- The durability and the solidity of the structure

The wheelchairs consists of the following elements:

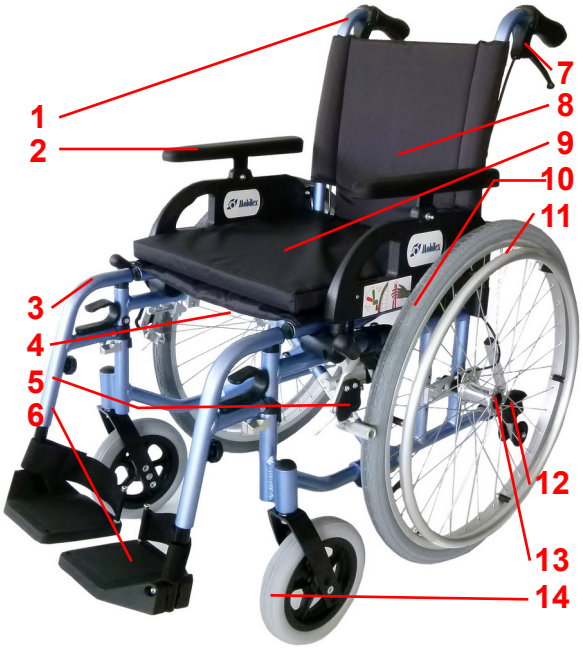
1. Push handles – enabling the attendant maneuvering the wheelchair.
2. Armrest – adjustable in height. It is used as a comfortable support of forearms. It also facilitates getting in and out of the wheelchair.
3. Footrest – The removable footrest is a safe and comfortable support for legs. The footrest is adjustable in height.
4. Cross (under the seat) – thanks to it the wheelchair can be easily fold for transportation and storage.
5. Brake – Pushing of brake lever towards front causes their pulling and protects the wheelchair against unforeseen movement and the rolling down the slope.
6. Footplate – User's feet rest on it. It is possible to lift up the footrest, which can ease getting in and out of the wheelchair.
7. Handle brakes (optional) – They enable controlling the wheelchair by an attendant.
8. Backrest – made of smooth incombustible material, enabling also easy upholstery cleaning.
9. Seat – Made in similar way to the backrest (with an optional seat cushion).
10. Main wheels – Majority of the user's weight rests on these 24 inches pneumatic/PU puncture proof wheels.
11. Push rim – Thanks to push rims fastened to main wheels self propelled drive of the wheelchair is possible.
12. Anti-tip (optional)– Prevents wheelchair from falling backwards.
13. Quick-release axle – With this function wheels can be easily taken off, which facilitates the transportation of the wheelchair.
14. Front wheels – 8 inches turning wheels are installed on rotating forks, which makes maneuvering of the wheelchair possible.

WARNING! Always pull both brakes, if you do not move on the wheelchair.
WARNING! Do not put fingers between the wheel and the push rims. It threatens the injury of the palm.

Function description

Footrests

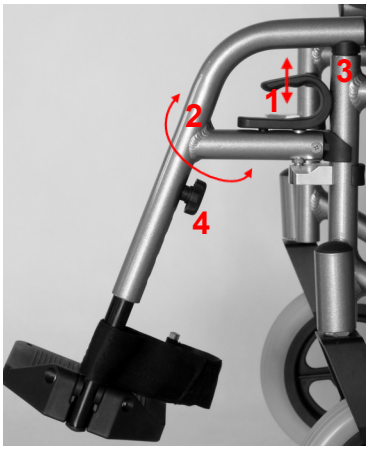Standard footrests can be disassembled or pushed away on side.

To rotate and disassemble the footrest, first unblock it by pulling up the safety-lock lever(1), then turn the footrest outside (2) and pull it upwards (3). The footrest assembling should be performed in reverse order by inserting it in base of frame (3) and turning it frontwards. The safety-lock will block automatically.


The wheelchair is equipped with height-adjustable footrests. To adjust their length one should loosen a knob (4), adjust the length of the footrest by sliding the lower tube up or down, and then retighten the knob.

During using the wheelchair outside, the footrests must be placed on height of minimum 4-5cm from the base. In addition, the footrests posses the ability of moving the footplates upwards, which facilitates getting in and out of the wheelchair and also disassembling the footrests.

WARNING! You should never climb the footrests, it can cause falling over the wheelchair. You should never raise the wheelchair using the footrests.

Mobilex A/S
Grønlandsvej 5
DK-8660 Skanderborg
www.mobilex.dk

Armrests

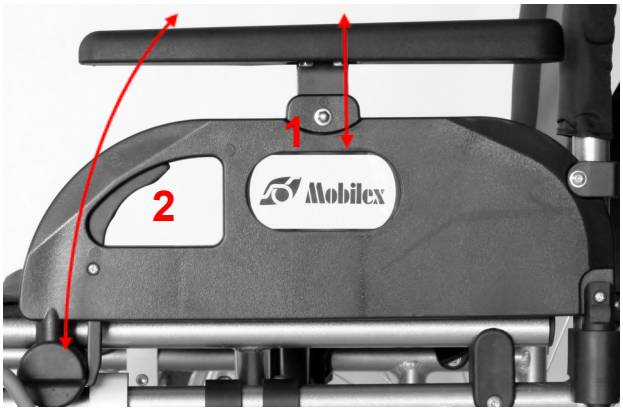The armrests can be adjusted in height, as well as be tilted backwards, behind the backrest, to facilitate the transfer of a passenger from or onto the wheelchair.


To adjust the height of the armrest, unscrew and remove the armrest bracket mounting screw (1). Adjust the armrest to the desired height, then insert and tighten the mounting screw back into the armrest bracket.

To tilt the armrest backwards, one should press the lever (2) and lift the armrest up.

To reattach the armrest, simply push it back down until the safety-lock will block itself automatically.

WARNING! It is prohibited to use the armrests as wheelchair carrying handles.

Parking Brakes

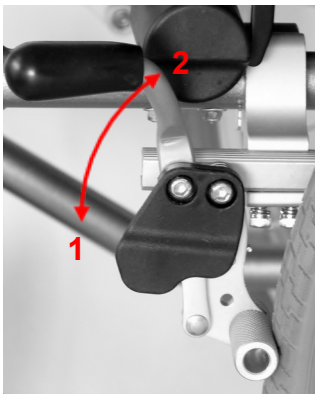
The wheelchair is equipped with parking brakes on both rear wheels to prevent rolling away.

To apply the brake push the lever forward (1). To release the brake pull the lever toward the wheel (2)

Always apply the brakes on both sides of the wheelchair.

WARNING! The parking brakes are not designed as driving brakes and only serve to lock the wheelchair (park position). They may not be used for braking while driving!

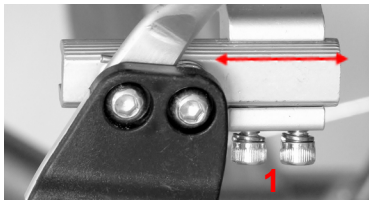The brakes can be adjusted for different types of drive wheels.

To set the correct distance between brake and tire:

Using a cross head screwdriver loosen the 2 screws (1) that secure the brake to the underside of the guide bar.
Slide the brake backwards or forwards to the correct position (approximately 10mm between the tire and the released brake shaft).
Then tighten the 2 screws.

Handle Brakes *(optional)*

The attendant brake levers are optional and may only be installed by the manufacturer or qualified service personnel. They assist the attendant in maneuvering the wheelchair user while driving.

The set consists of two brake levers (1) which control the wheelchair brakes and two steel cables (2) which connect the levers to the brakes. Brake levers are mounted under the push handles of the wheelchair.

To slow down the wheelchair, pull both brake levers (1) up evenly to control the braking force.

For full braking, pull the levers all the way up. To disengage the brakes, simply release the brake levers.

You can use the brakes to maneuver the wheelchair by braking the wheel on the side you want to turn to.

NOTE: The brake levers at the push handles may only be operated by persons assisting the user, standing behind the wheelchair. These levers are not intended to be operated by persons sitting in the wheelchair.


NOTE: The brake levers on the push handles are for maneuvering purposes and are not to be used as parking brakes.

NOTE: Never leave the wheelchair unattended without applying the parking brakes.

NOTE: Always apply the parking brakes when the wheelchair is not in use.

Mobilex A/S
Grønlandsvej 5
DK-8660 Skanderborg
www.mobilex.dk

24” rear wheels

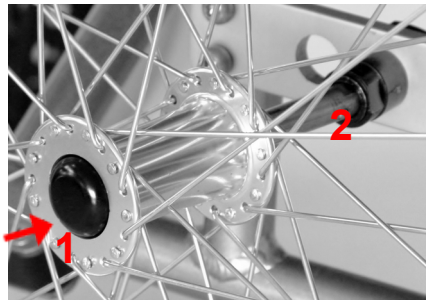The wheelchair is equipped with 24” main wheels with pneumatic or puncture-resistant polyurethane tyres. The wheelchair is delivered with the standard push rims. The material and the distance between the push rim and the wheel can influence on the grip strength of the wheelchair’s user.

One should be careful and not allow the fingers to get between spokes during driving in narrow passages. To avoid such risk it is recommended to use a spoke protection.

One should regularly control the pressure in tyres, which ought to be in range of 2 – 2,2 atm depending on the weight of the user.

The wheels on “Flipper” chair can additionally be easily removed and fitted in by pressing the quick-release button placed in the center of the wheel (1).

To install the wheel:

Place the end of the wheel axle back into the socket (2), press and hold the button (1), push the axle as far as it will go into the socket and release the button. Check if the wheel is fastened by pulling it outwards several times.


The wheel balance on the “Flipper” can be additionally adjusted by moving the wheel axle. In this case it is necessary to adjust also the position of the turning wheels.

WARNING! We recommend such adjustment to be done by the qualified staff. An inappropriate adjustment can cause deterioration of the safety usage of the wheelchair and decrease its durability.


Safety belts

There is possibility to attach the hip safety belts. They are not among the standard equipment.

How to fold and unfold the wheelchair

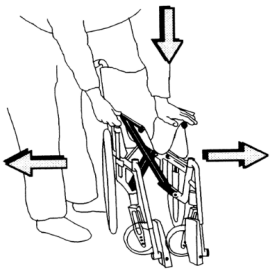Placing one hand on the seat upholstery press one of the tubes, which support the seat upholstery (in accordance with the instructions on the drawing). One should be careful and not put the hand between tubes and the wheelchair frame during unfolding it.

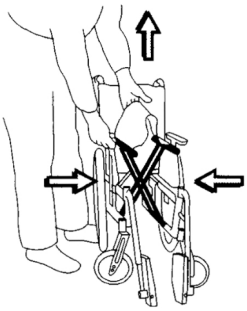Tip pedals

The tip pedal on wheelchairs makes it easier for the attendant to avoid obstacles by tilting the wheelchair backwards (see photo left). The tip pedal is built into the frame (see photo right).
NOTE! DO NOT tilt the wheelchair with the handles, but use the tip-pedal!

Advanced settings and adjustments

WARNING! We recommend that all advanced adjustments are only carried out by qualified personnel. An unsuitable adjustment may affect safety when using the wheelchair and reduce the service life of the wheelchair.


WARNING! The correct seat depth should always be professionally determined by a therapist or specialist!

Frame adjustments - changing the position of the rear and front wheels

The construction of the wheelchair makes it possible to change the height of the seat and also to shift the center of gravity of the wheelchair.

- Adjusting the seat height is done by changing the mounting position of the rear wheels (up/down) and changing the mounting position or size of the front wheels.
- Adjusting the balance of the wheelchair is done by changing the mounting position of the rear wheels (front/rear).

Adjusting wheel position (vertically and horizontally)
If necessary, the seat height of the wheelchair can be adjusted by raising or lowering the axles of the rear wheels. In this case, the position of the front wheels must also be adjusted.

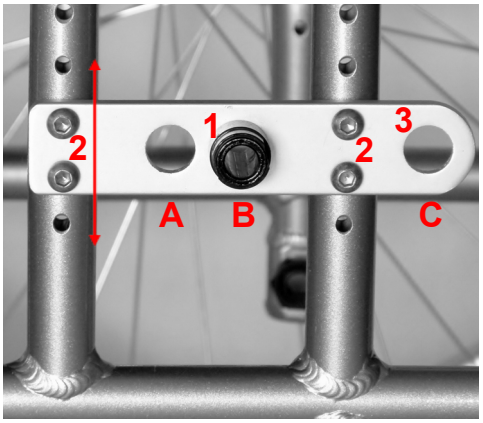It is possible to adjust the rear wheels up or down, forwards or backwards to suit the tyre size and the personal preferences of the user. The height adjustment of the wheel axles allows the adjustment of the seat height. The horizontal alignment (front / back) allows to influence the stability and manoeuvrability of the chair. The most common and neutral position is "B". If it is necessary to increase the manoeuvrability of the wheelchair, position "A" can be used. With the wheels in position "A" it is easier to manoeuvre the wheelchair but the stability to the rear is significantly reduced.
If the wheels are set to position "A" the chair should always be equipped with an activated anti-tip. Persons with amputated leg(s) should never use the wheelchair with the wheels in position "A" - the risk of tipping is too high. Leg amputees should set position "C" as it provides greater stability than the other positions. Even when using this setting, the wheelchair should be equipped with anti-tippers (accessories). As the manoeuvrability of the wheelchair is made more difficult with the axles in position "C", this setting is only recommended for leg amputees.

WARNING! The correct position of the wheels should only be determined by professional therapists or specialist companies.

There are 4 different height settings and 3 different horizontal settings possible.

To adjust the axles horizontally (front or back) proceed as follows:

1. Remove the wheels with the quick-release axle.
2. Remove the wheel mounting bracket (1) and place it in the desired position (A, B, C)
3. Repeat the process with the opposite axle.

To adjust the axle height, proceed as follows:

1. Remove the wheels with the quick-release axle.
2. Unscrew the wheel mounting plate by loosening the 4 screws (2).
3. Move the plate up or down to your desired height dictated by screw slots.
4. Fix the axle mounting plate again with the 4 screws.
5. Repeat the procedure with the opposite mounting plate.

8" front wheels

The wheelchair is equipped with PU wheels 200x50 mm as standard. We recommend to stay with this type of tyres, but it is also possible to use pneumatic tyres. It is possible to mount wheels with different sizes. This might be necessary after adjusting the rear axles. If the vertical position of the rear wheels is changed, it is always necessary to adjust the position of the front wheels accordingly. The vertical axis of rotation of the fork should always be set at a 90° angle to the ground (+/- 2°).

To make the adjustment, move the front wheel axle to another hole by loosening the bolt (1), moving the wheel up or down to the desired height and tightening it again in the selected hole. There are three holes available. In some cases it may be necessary to replace the standard front wheels (8") with a wheel of a different size (6" or 7").


Seat cushion (*optional*)

The wheelchair can optionally be equipped with a 3 cm thick seat cushion.

- Always remove the seat cushion from the seat before folding the wheelchair frame.
- After unfolding the frame, put the cushion back on the seat.
- Always position the cushion on the seat so that the cushion velcro tabs (1) are parallel and tightly connected with the velcro tabs on the seat.
- The zipper (2) allows you to easily remove the upholstery from the cushion in order to wash it.

Anti-tip support *(optional)*
This component increases rear stability and safety. It is especially recommended for wheelchairs of which center of gravity has been shifted rearwards. It should be installed on all wheelchairs used by leg amputees and is recommended for wheelchairs with the rear wheels mounted in the forward position.

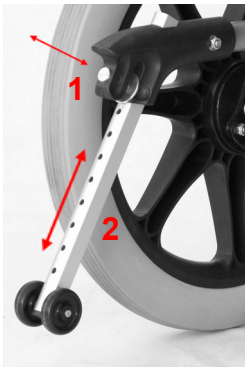To activate the support, pull and hold the locking pin backwards (1), then move the support (2) to a desired position. Release the locking pin (1) so that it locks the support in this position.

To deactivate the support, simply move the support all the way back after pulling the pin and let go once adjusted.

The safe distance between the support and the ground should not exceed 4 cm.

WARNING! Proper installation of the anti-tip support should only be performed by qualified personnel.

Assembly
The only elements removed from the device for delivery are the footrests. Insert them according to the **Footrests** section in **Function description** above.


Accessories
As an option, the wheelchairs can be equipped with the following accessories (not part of the standard equipment)

- ✔ Anti-tip
- ✔ Seat cushion
- ✔ Safety belt
- ✔ Handbrakes (on push handles)

4. Technical data

Flipper with / without optional equipment	**Item no. 271740 / 271940**	**Item no. 271744 / 271944**	**Item no. 271748 / 271948**	**Item no. 271751 / 271951**
Seat-width	40 cm	44 cm	48 cm	51 cm
Seat-depth adjustable	44 cm			
*Seat-height min./max.**	48-52,2 cm			
Total width (depending on choice)	61,5 cm	65,5 cm	69,5 cm	72,5 cm
Weight (depending on size and accessories)	up to 16-18 kg			
Max. user weight	136 kg			

Product label
A product label is placed on the front of the frame.

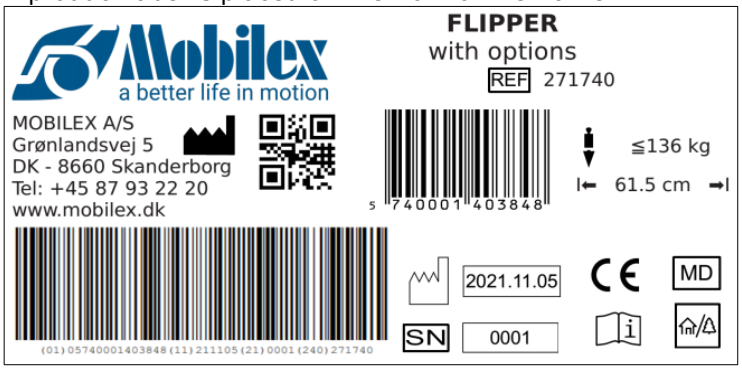Example of product label

	QR-code link to homepage		Manufacturer
CE	CE marked product	REF	Item Number
	Observe the user manual		Max load
MD	Medical Device		Manufacturing date
	To be used indoor and outdoor	SN	Serial no.
x cm	Max width		

Meaning of the icons


5. Care and maintenance

- ✔ The frame ought to be cleaned with use of damp cloth alternatively with addition of a soft detergent.
- ✔ **CAUTION!** The frame should be regularly controlled, to find on time the damages of varnish, which my cause corrosion.
- ✔ The seat and the backrest are made of a very durable strengthened nylon. The seat upholstery can be easily removed by slackening the screws placed on the frame tubes. The backrest upholstery can be removed by slackening the screws placed on the backrest tubes.
- ✔ The seat and the backrest upholstery can be washed with use of sponge and gentle washing soap.

Recommended service life: Normally 5 years depending on use.

Storage: We recommend storage in a dry place in temperatures above 0 °C.

Transportation
Place the wheelchair in the car boot. If, for some reason, it is not possible, put the folded wheelchair between the front and the rear seats inside the car ensuring that the wheelchair is in a sure and stable position. If it is not possible the wheelchair should be fasten with the safety belt.

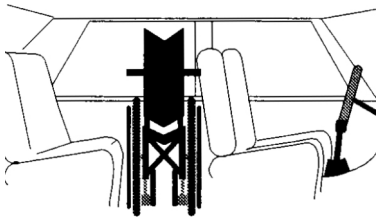WARNING! The wheelchair is not the car seat. During traveling by vehicles seating on the wheelchair is not allowed. The chair should be folded and safely immobilized.

Reuse
The product is suitable for reuse. Before a new user takes over the product, it must be checked according to point 5.

Parts, which should be regularly inspected

Part name	Control type	Control frequency
Tyres	Tyre pressure, condition of tread and tyres. One should take care, to keep the pressure in both tyres on the same level. Harder tyres provide better maneuverability of the wheelchair and larger easiness of driving, but they reduce the comfort on uneven surface.	At least once a week
Spokes	Loose spokes can cause deformation of the rear wheels. Contact your point of sale or local bicycle service to manage the problem of loose spokes.	If the problem occurs
Wheel axles	Remove the hair or accumulated dirt	If necessary
Push Rims	Excessively scratched push rims ought to be exchanged because they can wound the user’s hands during driving the wheelchair	If necessary
Brakes	The braking force depends on the tyre pressure. The efficiency of brake operation can be also affected with dirt accumulated on tyres. Keep the brakes clean wiping them with a damp cloth to remove dirt, and also lubricating the funnel of screws, on which the brake levers turn.	Control the correctness of brake operation at least once a month
Frame	Maintain the wheelchair in the highest cleanness for better comfort of the user.	At least once a month
Turning wheels	The area between the fork and the front wheel should be kept clean, because of dirt accumulating there, which can cause faster wear of turning wheel bearings. To do so, one should disassemble the front wheel by unscrewing it from the fork, remove all dirt, and then preserve the metal elements of the wheel (i.e. technical grease).	The maintenance ought to be made once a month or often depending on the soiling degree.
Detachable elements	Check condition of detachable elements of the wheelchair, if screws are loosen they should be tighten.	In case of intensive exploitation of the wheelchair the inspection should be carry out once a month.
Quick-release axles of the rear wheels	One should clean and then lubricate lightly the quick-release axles of the rear wheels, with a damp cloth the slots of fixing the quick-release axles ought to be cleaned up from all dirt. The wheel should be laid on the base, resting on the push rim (never opposite way, because the quick-release axle can be soil, which can cause its seizing in the main axle opening. In case of solid particles occurring on the surface of the quick-release axle of the wheel, using is not allowed.	Lubricating should be carried out periodically.

6. Most common problems and solutions

If you notice any irregularities in the wheelchair working you cannot use it. You should contact the local point of sale or service department of the wheelchair supplier. The manufacturer does not guarantee correctness of the wheelchair operation, if the parts applied in it are not original.

Problems	*Solutions*
The wheelchair seems to tilt to one side	• inflate tyres • check fastening of the front turning wheels • check displacement of the wheelchair mass • one of the rear tyres might be inflated more then the other one
The wheelchair is hard to push	• inflate tyres • check fastening of the front turning wheels • remove dirt or entangled hair from the front turning wheels • to large load applied on the front turning wheels: move the center of gravity
The wheelchair is hard to turn	• inflate tyres • ensure, that bolts fastening the front turning wheels are not tightened to strong • remove dirt or entangled hair from the front turning wheels
Brakes are not working correctly	• inflate tyres • adjust the distance between the brake and tyres
It is hard to fold and unfold the wheelchair	• the upholstery fitted to tight • clean up or besmear the cross under the seat
The wheelchair is not stable	• inflate tyres • ensure, that all screws/bolts and nuts are tightened • ensure, the wheelchair is properly adapted to the passenger

For realization of basic service actions the set of socket wrenches, or adjustable spanner, flat and cross screwdriver will be sufficient.

7. Disposal of the product

The product can not be disposed of with household waste but must be brought to the local recycling center.

8. Guarantee terms

Mobilex A/S offers 2 years of warranty for damages caused by errors in production or material errors. The warranty is not valid by not-intended repair or usage of the product. Parts subject to normal wear are not covered by the warranty, unless wear is caused by a manufacturing fault. If product changes are made without our written acceptance the CE-marking and product warranty will not be valid. In the event of product damages covered by this warranty, please notify your dealer or Mobilex A/S directly. The warranty does not cover transport costs and does not include compensation for personal injuries or non-use of the product during repair. The warranty does not cover damages due to the user manual not being followed.

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Mobilex A/S
Grønlandsvej 5
DK-8660 Skanderborg
www.mobilex.dk

Mobilex A/S, Grønlandsvej 5, DK-8660 Skanderborg, Denmark, SRN: DK-MF-000021885

ENGLISH
Declaration of Conformity

Hereby declare under our sole responsibility as a legal manufacturer that the product specified on the product list below, meet the essential health and safety requirements and is in conformance with the provisions of the Regulation (EU) 2017/745 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL on medical devices.

The product are classified as Class I, medical device. The classification is based on the requirements of Rule 1 of annex VIII, of the Regulation (EU) 2017/745.

The CE marking has been affixed on the product according to Annex V of the Regulation (EU) 2017/745.

Intended purpose
Indicated for people with mobility limitations who have problems with walking, standing or sitting without stable support

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